

en*Puls* 2.0  
the next generation



en*Puls* 2.0 shown with optional aluminum transport case



en*Puls* 2.0 shown with optional *move* cart and hand piece holder

**Zimmer**  
MedizinSystems  
3 Goodyear, Suite B  
Irvine, CA 92618  
Phone: (800) 327-3576  
info@zimmerusa.com  
www.zimmerusa.com  
shop.zimmerusa.com

en*Puls* 2.0 Specifications

Technology	Compressor-free ballistic Radial Pulse Therapy system with electromagnetic generator for projectile acceleration
Power Levels	60-185 mJ freely selectable in 10 mJ increments
Modes	Frequencies from 1 to 22 Hz 3 Burst modes (4 / 8 / 12 pulses)
Programs	7 preset programs, customizable
Controls	7" color touchscreen for all software operations Rotary control for energy level and for frequency Hand piece operation with multi-directional foot switch
Memory / Updates	SD card for favorites and customized programs (up to 120 entries) Firmware updates
Radial Pulse Hand Piece	Ergonomic, with anodized aluminum housing and cooling fan
Generator	Electromagnetic, integrated
Dimensions	9.1" L, 2" diameter (max.)
Weight	ca. 1.87 lbs. (with cable)
Warranty	2,000,000 pulses
Applicator heads	Diameters of 6 / 15 (2 pieces) / 25 mm Changeable without tools
Power Consumption	100 - 240 V / 50/60 Hz, max. 250 VA
Dimensions	12.7 x 9.3 x 5.1 inches (L/W/H)
Weight	6 lbs. (without hand piece)
Package	Foot switch Mains cable Applicator head (6 / 15 (2 pieces) / 25 mm) 10 blue silicone caps 1 hand piece 1 bottle en <i>Puls</i> lotion
Options	Aluminum transport case <i>move</i> cart Holder for 1 hand piece

USA V2 | In the interest of product improvement, the specifications and/or appearance of products are subject to change without notice. | Zimmer MedizinSystems is in no way affiliated with Zimmer Biomet™

**Zimmer**  
MedizinSystems

en*Puls* 2.0  
Radial Pulse Therapy

Zimmer MedizinSystems is in no way affiliated with Zimmer Biomet™



**enPuls 2.0**

the next generation

enPuls & move



enPuls 2.0 shown with optional move cart and hand piece holder

## Radial Pulse Therapy

An Effective Alternative For Treating Chronic Soft Tissue Pain

### Principle of Radial Pulse Therapy

The hand piece contains a projectile that is accelerated through the electromagnetic transfer of kinetic energy. This kinetic energy is then transformed into impact energy in the applicator head. The impact energy delivered from the applicator head results in radial pulses in the target tissue. Radial Pulse Therapy provides an excellent alternative to treating muscle tissue. The energy treats muscle pain and helps to increase circulation in the target tissue.

### Mechanical energy

The enPuls device delivers energy into the target tissue where traditional soft tissue therapy techniques cannot reach. This is all done with little or no stress on the therapist's hands. A heavy duty, ergonomic aluminum hand piece delivers optimal energy to the patient while providing optimal protection against vibrations to the therapist's hands.

### Time savings and efficiency with better outcomes

Adding the enPuls Radial Pulse Therapy unit to your facility will decrease the treatment time for chronic soft tissue pain patients and increase your efficiency. This results in better outcomes for you and for your patients.

## Stimulating and gentle on the patient...



Exceptionally Powerful Technology

### enPuls Softshot technology

enPuls creates unique radial pulses by using a projectile 5 times heavier than standard air compressor systems and a reduced final velocity upon contact with the applicator. Though the force of the accelerator and the mass of the projectile are efficient down to 35 mm depth, they are more comfortable for the patient than comparable systems.

### Increased patient compliance

Less pain, same energy – that's enPuls softshot technology. Energy is maintained by a longer pulse rise time. This eliminates the need for unpleasant high peak pressure amplitudes. The lower peak and longer rise time minimize the stinging and painful sensation that is normally associated with this type of therapy. Increased patient compliance results in easier treatments at higher energy levels, leading to successful therapy.

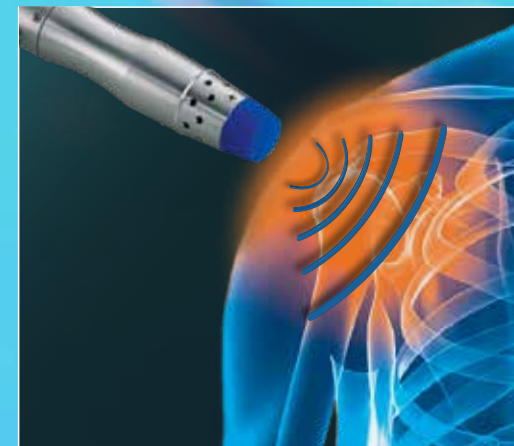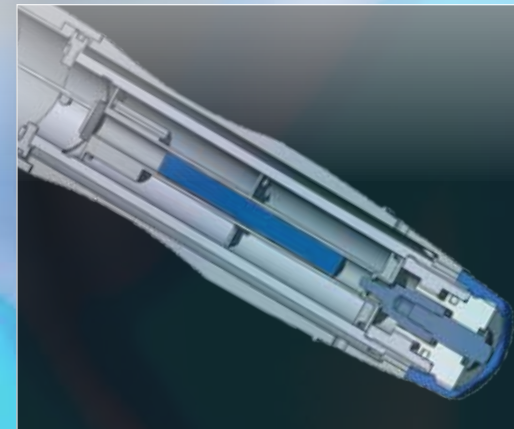

**Zimmer**  
MedizinSystems



# enPuls 2.0

the next generation

## Details

- 1 hand piece
- intuitive touch-screen operation
- energy levels from 60 to 185 mJ (comparable to 1-5 bar)
- pulse frequencies from 1 to 22 Hz
- burst mode for trigger point treatment
- positive/negative pulse counter with value preset
- SD card for memory, updates and service
- efficient energy transfer
- lower maintenance, longer lifetime
- unique, comfortable pulses (softshot)

## ... efficient and easy to use

### Ergonomic

The high mass applicator reduces strain on clinicians. The ergonomic aluminum hand piece with high mass and ample inertia provides appropriate energy transfer to the patient. It also affords optimal protection against vibration for the healthcare professional. The design allows for single handed use allowing the second hand to be free to adjust treatment parameters in session.



### Safe

Pulses are triggered using a foot pedal, ensuring targeted treatment and allowing the clinician to concentrate on hand positioning for the treatment.



### Flexible

Applicator heads of varying sizes and shapes can be quickly changed out without the use of tools, allowing treatments to be adjusted as needed.



### Mobile

An optional aluminum transport case containing all equipment and accessories offers mobility and convenience for treatment any time and anywhere.



## Revolutionary interface design

The new enPuls touch screen interface is truly innovative, intelligently designed... and even fun to use.



## Larger display surface

The 7" color touchscreen display is easy to read and to navigate with your fingertips.

### 2 Million Pulse Warranty

A powerful, mechanical energy generator with high percussion, digital control, and low mechanical wear provides for an extended lifetime of reliability and performance. We guarantee our hand piece for 2 million pulses.

### No Air Compressor

The enPuls is completely portable for use between treatment rooms or on the road. There are no messy air compressors to drag, listen to or maintain.



**enPuls 2.0**

the next generation

## Successful Therapy

Radial Pulse Therapy: easy to operate with quick results

*The enPuls becomes an extension of your hands, delivering mechanical stimulation to the target tissue.*

Radial Pulse Therapy provides an excellent alternative for treating muscle tissue.

The energy treats muscle pain and help increase the circulation in the target tissue.

Easy-to-use and requires minimal preparation, enabling the practitioner to deliver efficient and effective treatments.

On average, Radial Pulse Therapy is performed once or twice a week for a total of 5-6 sessions.

Most treatments are administered at 2,500 - 3,000 pulses and take less than 10 minutes per session.



**Zimmer**  
MedizinSystems

## Treat soft tissue pain in your office with this cost-effective, non-surgical alternative

It's had a tremendous effect on our practice, enabling patients to get better sooner.

*– D. Lipitz, Physical Therapist*

Pulsed wave therapy is a tremendous asset to my practice. It is the most effective treatment modality that I use to treat patients suffering from chronic foot pain unrelieved by traditional conservative treatments.

*– M. Cunha, Doctor of Podiatric Medicine*

Nothing I have used has been as significant as Radial Pulse Therapy. It saves my hands.

*– T. Rodman, Chiropractic Physician*

After two back surgeries and one spinal cord pain stimulator did not relieve my pain, the doctors said I would just have to live with it. After RPT I feel I have recovered at least 80% of my quality of life.

*– J. Piescharon, 73-year-old Patient*

